

Epic Radiant

Recent Studies Report with Result History

Custom Reports

April 1, 2025

Epic Build Description

For each exam on the Recent Studies Report, we display Study Details, as shown on the right-side report.

In Epic, we changed the print groups that display on this report so that previous versions of the result report could be viewed.

The screenshot displays the 'Recent Studies Report' interface in Epic. The main table lists studies with columns for MRN, Procedure, and End Exam Time. A red box highlights the 'Recent Studies Report' tab. To the right, a 'Study Details' panel is open for the study 'MRI STEALTH WITHOUT CONTRAST (70551) [IMG3218]'. The panel shows the reading physician, reading date, and a narrative & impression section.

MRN	Procedure	End Exam Time
	MRV HEAD W/WO CONTRAST (70546)	Tue Mar 25, 2025 10:57 AM
	MRI STEALTH WITHOUT CONTRAST (...)	Fri Mar 28, 2025 2:48 PM
	MRI SPINE THORACIC W/O CONTRAS...	Fri Mar 28, 2025 12:02 PM
	MRI SPINE THORACIC W/WO CONTR...	Thu Mar 27, 2025 9:57 AM
	MRI SPINE THORACIC W/WO CONTR...	Thu Mar 27, 2025 9:03 AM
	MRI SPINE THORACIC W/WO CONTR...	Tue Mar 25, 2025 3:10 PM
Chuck	MRI SPINE LUMBAR W/O CONTRAST (...)	Wed Mar 26, 2025 3:27 PM
	MRI SPINE LUMBAR W/O CONTRAST (...)	Wed Mar 26, 2025 8:56 AM
Pat	MRI SPINE LUMBAR W/WO CONTRAS...	Fri Mar 28, 2025 10:36 AM
	MRI SPINE LUMBAR W/WO CONTRAS...	Tue Mar 25, 2025 3:10 PM
	MRI SPINE CERVICAL W/O CONTRAST...	Fri Mar 28, 2025 12:55 PM
	MRI SPINE CERVICAL W/O CONTRAST...	Fri Mar 28, 2025 12:02 PM
	MRI SPINE CERVICAL W/O CONTRAST...	Fri Mar 28, 2025 7:33 AM
Geoff	MRI SPINE CERVICAL W/O CONTRAST...	Wed Mar 26, 2025 9:23 PM
	MRI SPINE CERVICAL W/O CONTRAST...	Wed Mar 26, 2025 7:58 AM
	MRI SPINE CERVICAL W/WO CONTRA...	Fri Mar 28, 2025 9:28 AM
	MRI SPINE CERVICAL W/WO CONTRA...	Thu Mar 27, 2025 10:16 AM
	MRI SPINE CERVICAL W/WO CONTRA...	Thu Mar 27, 2025 9:04 AM
	MRI SPINE CERVICAL W/WO CONTRA...	Tue Mar 25, 2025 9:55 AM
	MRI ORBIT W/WO CONTRAST (70543)	Thu Mar 27, 2025 7:58 AM
	MRI BRAIN W/O AXIAL T2 ONLY (70551)	Fri Mar 28, 2025 10:02 AM
	MRI BRAIN W/WO CONTRAST (70553)	Fri Mar 28, 2025 4:25 PM
	MRI BRAIN W/WO CONTRAST (70553)	Fri Mar 28, 2025 1:20 PM
	MRI BRAIN W/WO CONTRAST (70553)	Fri Mar 28, 2025 11:45 AM
	MRI BRAIN W/WO CONTRAST (70553)	Fri Mar 28, 2025 10:18 AM
	MRI BRAIN W/WO CONTRAST (70553)	Fri Mar 28, 2025 8:25 AM
	MRI BRAIN W/WO CONTRAST (70553)	Thu Mar 27, 2025 10:16 AM
	MRI BRAIN W/WO CONTRAST (70553)	Thu Mar 27, 2025 9:03 AM

Study Details

MRI STEALTH WITHOUT CONTRAST (70551) [IMG3218]

Order #: 57 Accession #: RAD96

Reading Physician: MD
Reading Date: 3/28/2025
Result Priority: Routine

Narrative & Impression

Procedure: MRI STEALTH WITHOUT CONTRAST (70551)

Indication: Encounter for other preprocedural examination. preop mapping DBS implant 3T scanner - T1, T2, SWI, stealth VIZ DTI on 3T scanner

Technique: Axial SPGR volumetric and DTI images of the brain was performed without IV contrast. Exam was performed for treatment planning.

Comparison: None.

IMPRESSION


Findings/Impression:

- Exam was performed for pretreatment planning.
- No acute intracranial pathology.

Epic Build Description


Epic has a standard print group called “Results: History of Order Results” (LPG 52100). We made a copy of this print group and added it to our report. It displays the result history in reverse chronological order.

The link in the top right corner allows the option to toggle between result history view and result comparison view.

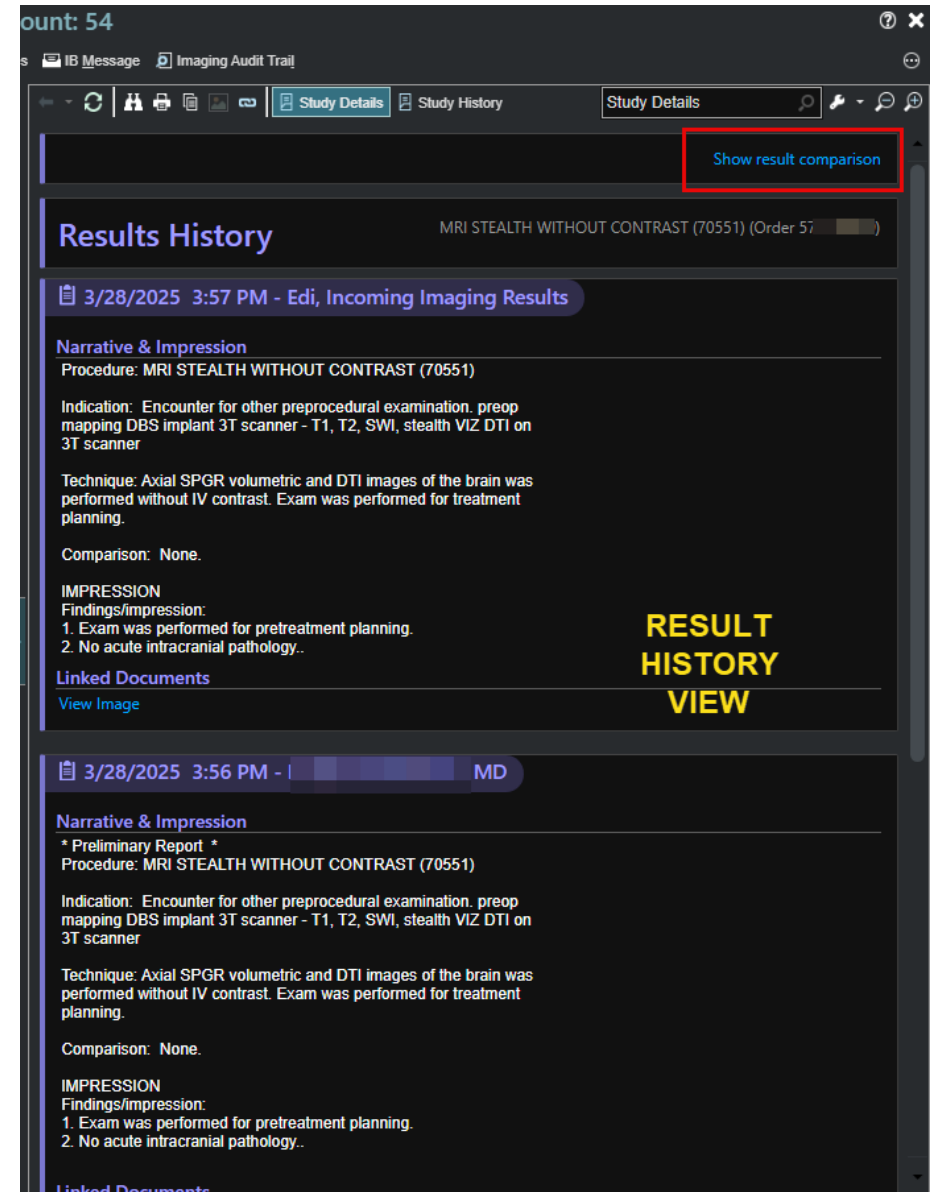


Epic Build Description

This is the result comparison view.
The print group can be customized to use either of these views as the default view.

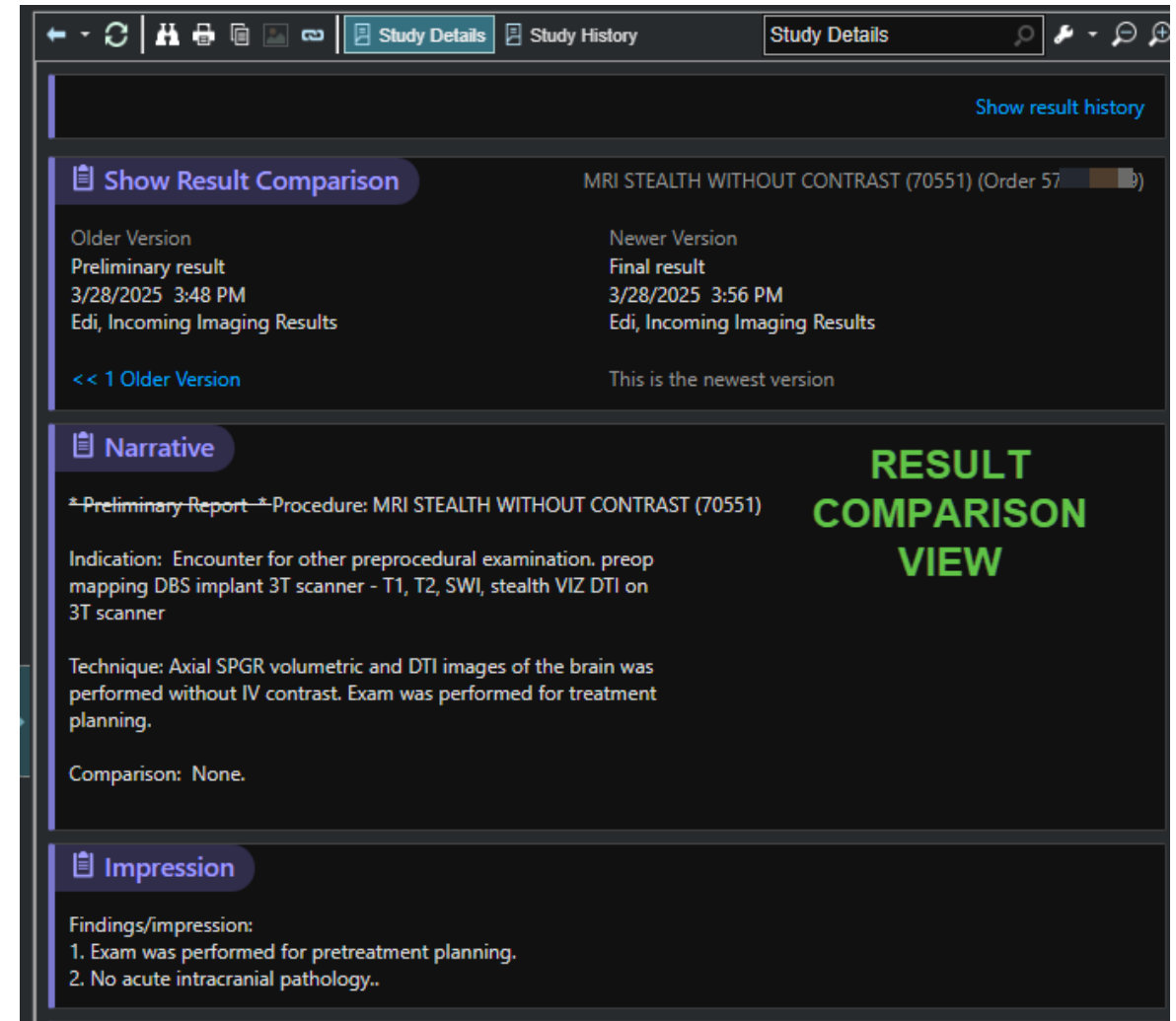

Epic Build Steps

1. Copy Print Group 52100 to create a custom print group.
 - a. In Epic Text, edit the new print group.
Change parameter 27. Default Display Mode to:
“0-List Mode” to display result history view as the default or
“1-Comparison Mode” to display result comparison view as the default

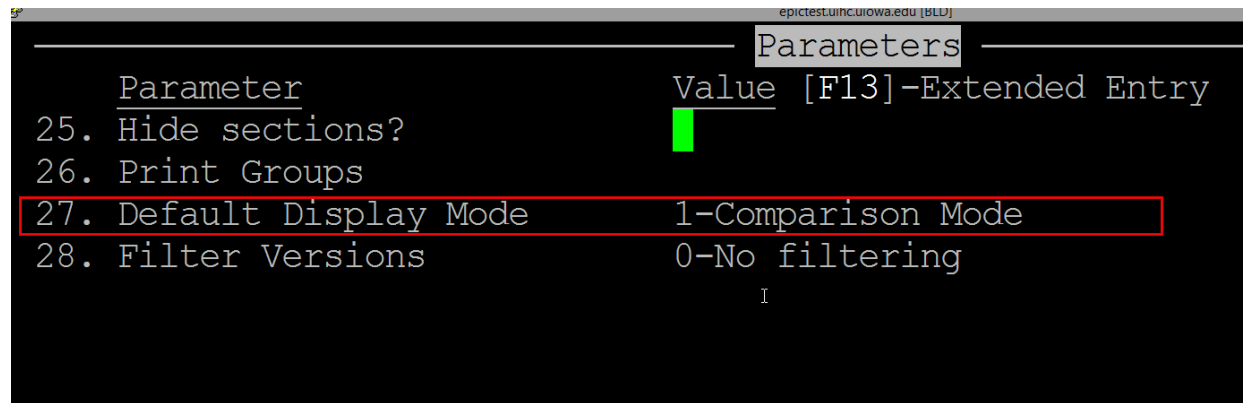

The screenshot shows a terminal window titled "epictest.uinc.iowa.edu [BLU]". Inside, a table titled "Parameters" is displayed. The table has two columns: "Parameter" and "Value [F13]-Extended Entry". The following parameters are listed:

Parameter	Value [F13]-Extended Entry
25. Hide sections?	
26. Print Groups	
27. Default Display Mode	1-Comparison Mode
28. Filter Versions	0-No filtering

The row for parameter 27 is highlighted with a red border. A green cursor is visible at the end of the value "1-Comparison Mode".

Epic Build Steps

2. Identify the currently used “Recent Studies report” on your profile record
 - a. If the current report on your profile record is a custom report, you can modify it. Proceed to step 3.
 - b. If the current report is an Epic-released report, make a copy of it. Then update the profile to use the new copy. LPR item 52100 for Recent Studies Report [7]
= <new copy of LRP>
3. Modify your report (LRP) to include the custom print group (LPG) you created in step 1.

IOWA

Questions?

→ uiowa.edu

Cass Garrett

Lead Application Developer-Radiant Team
Application Services
University of Iowa Health Care
cass-garrett@uiowa.edu